


Manawa^{FLEX}

URETERAL ACCESS SHEATH

SMOOTH INSERTION
FOR CONTINUOUS
URETERAL ACCESS

MINIMIZED TIME TO
PERFORM THE
PROCEDURE IN
ONE STEP

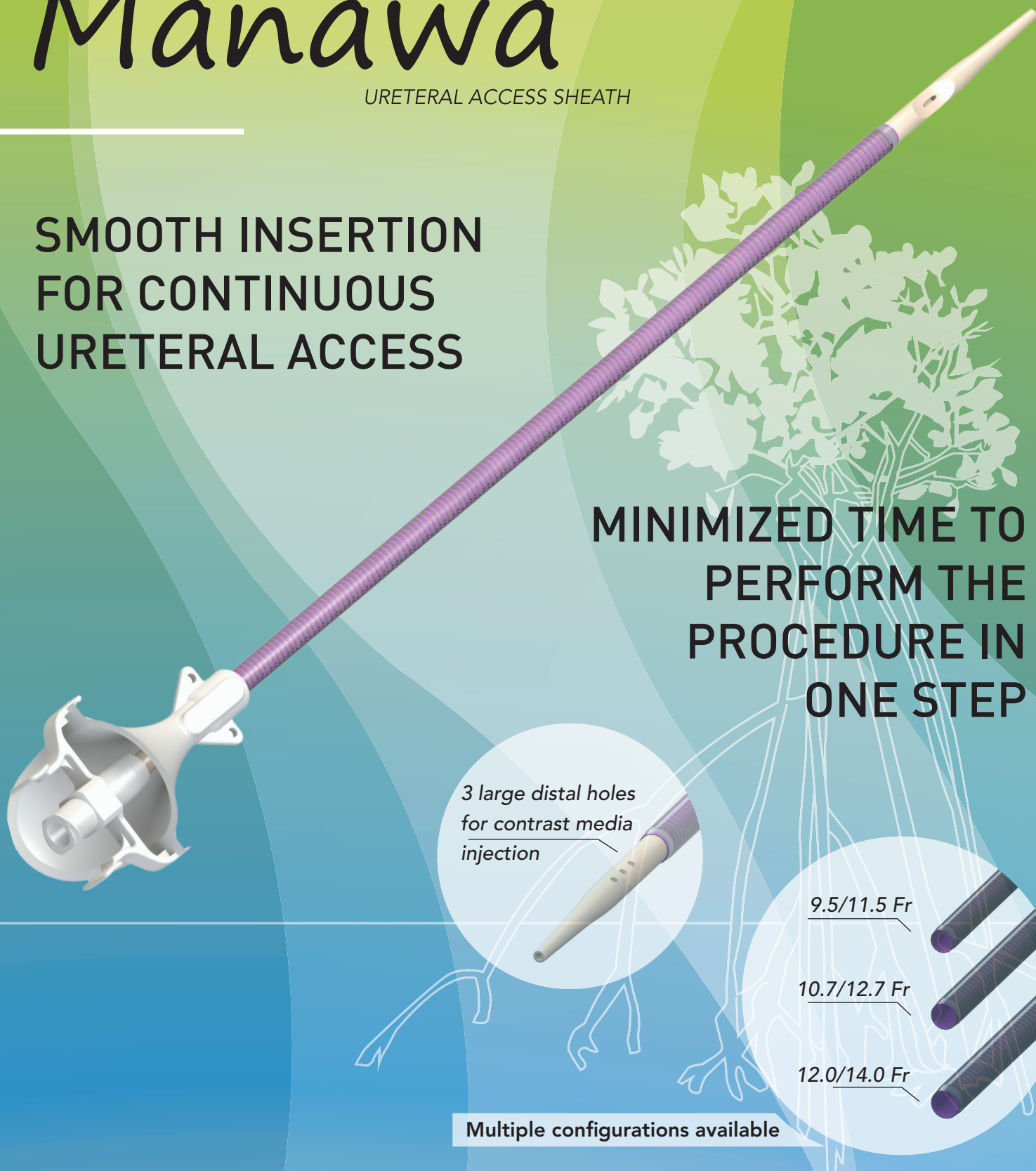


The Manawa Flex ureteral access sheath provides ureteral dilation with continuous working channel for the introduction of ureteroscopes and instruments during ureteroscopic procedure. It reduces the number of devices used: save one guidewire and one dual lumen catheter while allowing contrast medium injection during insertion.

Effective access for safer procedure

- Protect the ureter during repeated instrument exchange minimizing the risk of trauma
- The PTFE lined sheath protects delicate instruments and ureteroscopes from damage
- Faster procedure time : in one step the working guidewire becomes a safety guidewire
- Reduced number of devices: saves one guidewire and one dual lumen catheter
- Possibility to inject contrastmedia during the access sheath insertion

Smooth insertion and gentle dilation

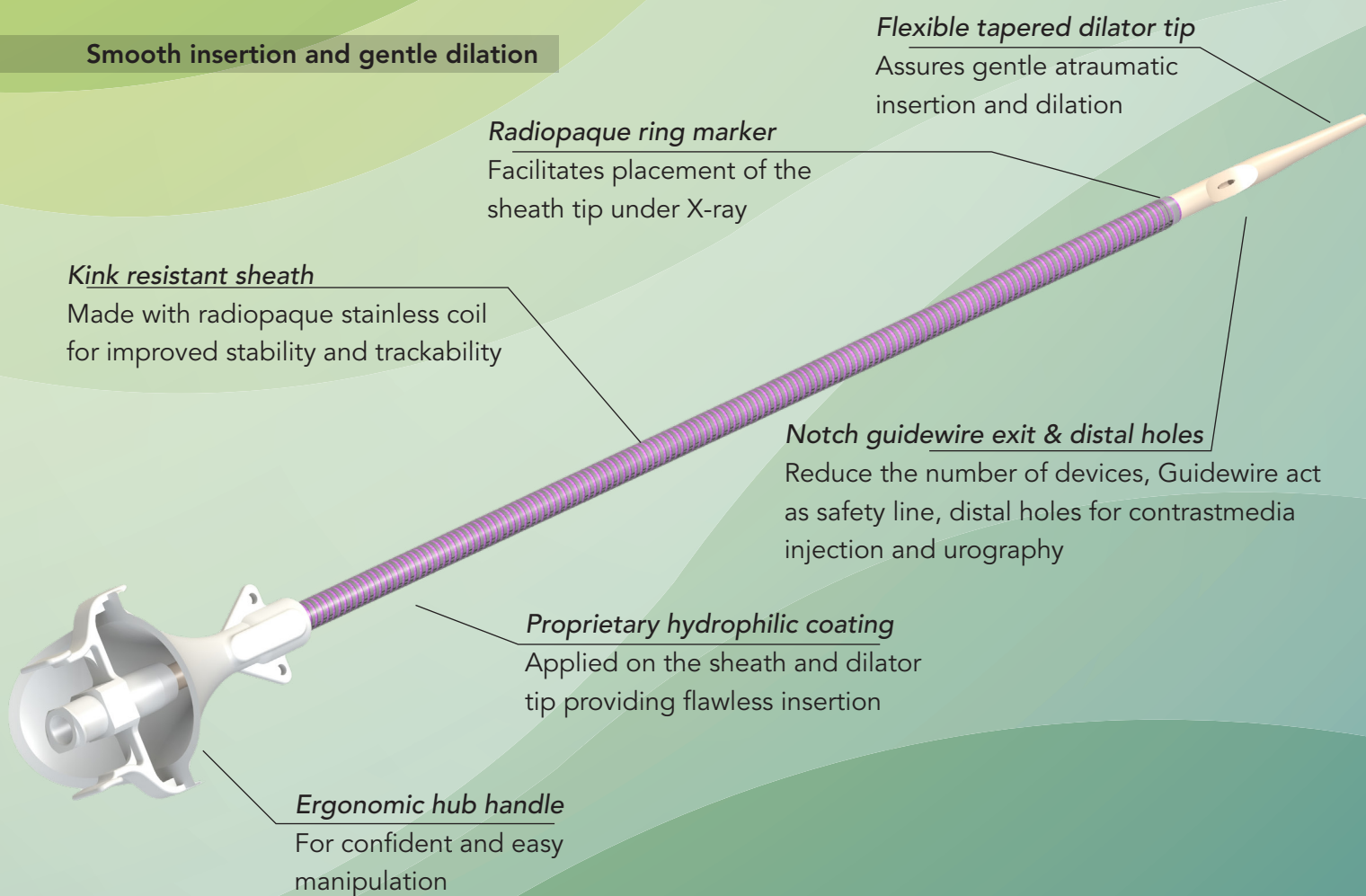

Order Number	Internal Diameter French Size	External Diameter French Size	Length
<i>Manawa Flex access sheath</i>			
UAS-0309520S	9.5 Fr	11.5 Fr	20 cm
UAS-0309528S	9.5 Fr	11.5 Fr	28 cm
UAS-0309535S	9.5 Fr	11.5 Fr	35 cm
UAS-0309545S	9.5 Fr	11.5 Fr	45 cm
UAS-0310720S	10.7 Fr	12.7 Fr	20cm
UAS-0310728S	10.7 Fr	12.7 Fr	28cm
UAS-0310735S	10.7 Fr	12.7 Fr	35cm
UAS-0310745S	10.7 Fr	12.7 Fr	45cm
UAS-0312020S	12.0 Fr	14.0 Fr	20 cm
UAS-0312028S	12.0 Fr	14.0 Fr	28 cm
UAS-0312035S	12.0 Fr	14.0 Fr	35 cm
UAS-0312045S	12.0 Fr	14.0 Fr	45 cm