


NEPHROSTOMY VACUUM

Manawa

NEPHROSTOMY ACCESS SHEATH WITH SUCTION

CONTINUOUS FLOW NEPHROSTOMY SHEATH FOR PERCUTANEOUS NEPHROLITHOTOMY (PCNL) PROCEDURE

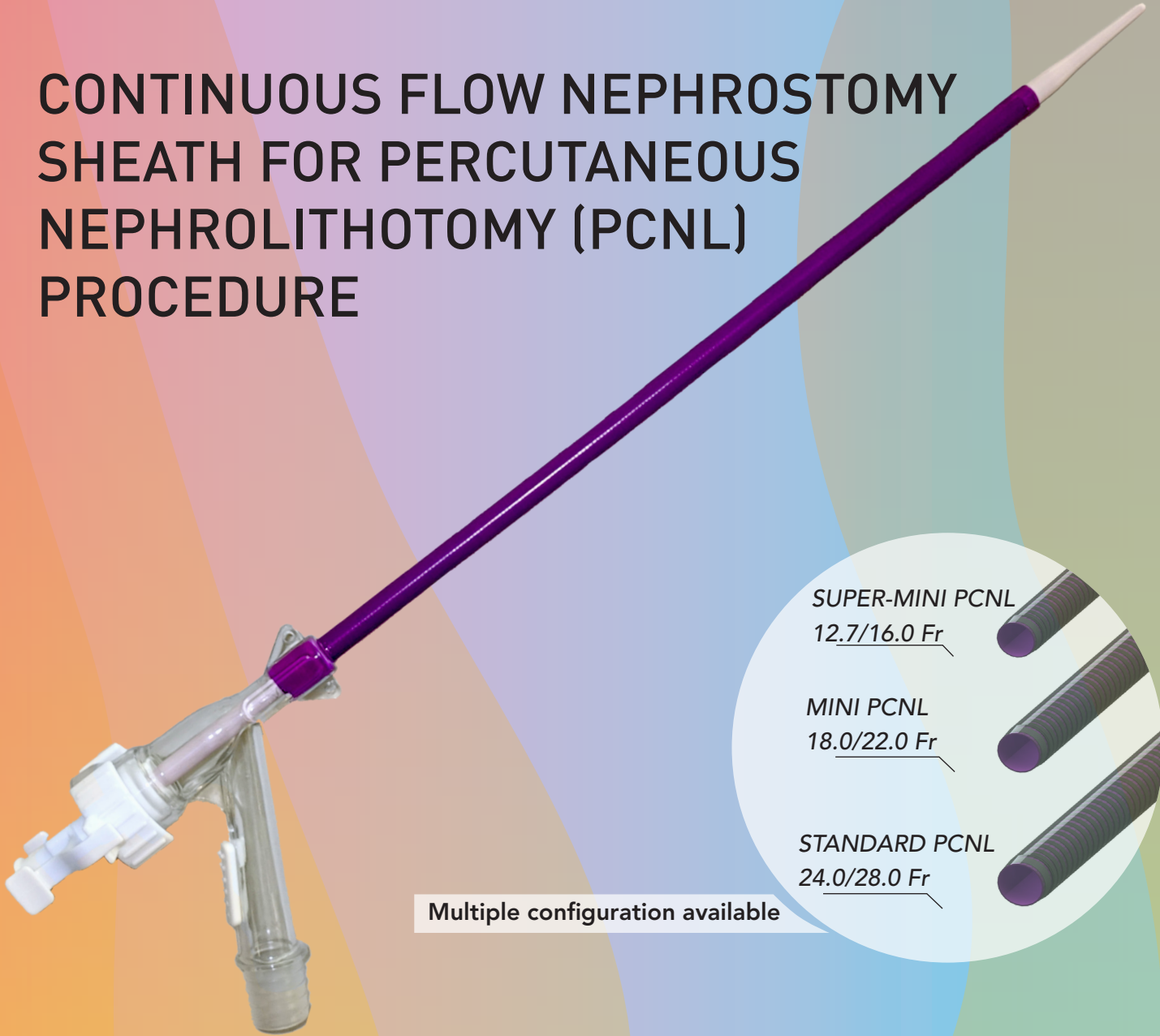

SUPER-MINI PCNL

12.7/16.0 Fr

MINI PCNL

18.0/22.0 Fr

STANDARD PCNL

24.0/28.0 Fr

Multiple configuration available

The Manawa Nephrostomy suction access sheath reduces intrarenal pressure, improve vision, remove stone fragments and dust while providing a safe, continuous working channel for the introduction of nephroscopes and instruments during PCNL procedure.

Effective access for safer procedure

- Reduces intra-luminal pressure thanks to a vortex created by the continuous irrigation and suction
- Superior pushability of the sheath to facilitate access the kidney calyx
- Unique funnel cone design handle facilitating scope entry, protecting your high-value optics
- Innovative valve design to ensure efficient suction control and negative pressure management
- Manual control of the suction intensity via pressure ventilation slider
- Hydrophilic coating on the dilator and sheath for smooth insertion
- Tapered dilator with rounded tip for atraumatic introduction
- Seamless transition between dilator and sheath to reduce risks of trauma

Order Number	Internal Diameter French Size	External Diameter French Size	Sheath Length	Hydrophilic Coating
Manawa nephrostomy vacuum access sheath				
SUPER-MINI PCNL				
UAS-0510713S	10.7 Fr	12.7 Fr	13 cm	dilator tip & sheath
UAS-0510715S	10.7 Fr	12.7 Fr	15 cm	dilator tip & sheath
UAS-0510717S	10.7 Fr	12.7 Fr	17 cm	dilator tip & sheath
UAS-0510721S	10.7 Fr	12.7 Fr	21 cm	dilator tip & sheath
UAS-0512013S	12.0 Fr	14.0 Fr	13 cm	dilator tip & sheath
UAS-0512015S	12.0 Fr	14.0 Fr	15 cm	dilator tip & sheath
UAS-0512017S	12.0 Fr	14.0 Fr	17 cm	dilator tip & sheath
UAS-0512021S	12.0 Fr	14.0 Fr	21 cm	dilator tip & sheath
UAS-0514013S	14.0 Fr	16.0 Fr	13 cm	dilator tip & sheath
UAS-0514015S	14.0 Fr	16.0 Fr	15 cm	dilator tip & sheath
UAS-0514017S	14.0 Fr	16.0 Fr	17 cm	dilator tip & sheath
UAS-0514021S	14.0 Fr	16.0 Fr	21 cm	dilator tip & sheath
MINI PCNL				
UAS-0516013S	16.0 Fr	18.0 Fr	13 cm	dilator tip & sheath
UAS-0516015S	16.0 Fr	18.0 Fr	15 cm	dilator tip & sheath
UAS-0516017S	16.0 Fr	18.0 Fr	17 cm	dilator tip & sheath
UAS-0516021S	16.0 Fr	18.0 Fr	21 cm	dilator tip & sheath
UAS-0518013S	18.0 Fr	20.0 Fr	13 cm	dilator tip & sheath
UAS-0518015S	18.0 Fr	20.0 Fr	15 cm	dilator tip & sheath
UAS-0518017S	18.0 Fr	20.0 Fr	17 cm	dilator tip & sheath
UAS-0518021S	18.0 Fr	20.0 Fr	21 cm	dilator tip & sheath
UAS-0520013S	20.0 Fr	22.0 Fr	13 cm	dilator tip & sheath
UAS-0520015S	20.0 Fr	22.0 Fr	15 cm	dilator tip & sheath
UAS-0520017S	20.0 Fr	22.0 Fr	17 cm	dilator tip & sheath
UAS-0520021S	20.0 Fr	22.0 Fr	21 cm	dilator tip & sheath
STANDARD PCNL				
UAS-0522013S	22.0 Fr	24.0 Fr	13 cm	dilator tip & sheath
UAS-0522015S	22.0 Fr	24.0 Fr	15 cm	dilator tip & sheath
UAS-0522017S	22.0 Fr	24.0 Fr	17 cm	dilator tip & sheath
UAS-0522021S	22.0 Fr	24.0 Fr	21 cm	dilator tip & sheath
UAS-0524013S	24.0 Fr	26.0 Fr	13 cm	dilator tip & sheath
UAS-0524015S	24.0 Fr	26.0 Fr	15 cm	dilator tip & sheath
UAS-0524017S	24.0 Fr	26.0 Fr	17 cm	dilator tip & sheath
UAS-0524021S	24.0 Fr	26.0 Fr	21 cm	dilator tip & sheath
UAS-0526013S	26.0 Fr	28.0 Fr	13 cm	dilator tip & sheath
UAS-0526015S	26.0 Fr	28.0 Fr	15 cm	dilator tip & sheath
UAS-0526017S	26.0 Fr	28.0 Fr	17 cm	dilator tip & sheath
UAS-0526021S	26.0 Fr	28.0 Fr	21 cm	dilator tip & sheath



Maeva Centre Building
Silicon Avenue, Ebene Business
Park, Reduit 72201
Mauritius

T: +230 466 3054